Features of New **Dr. SOS+ Kit**

01. Solve all cases with one kit

Capable of solving all cases such as fixture fracture, abutment fracture, screw fracture, 1.2 hex deformation, etc.

02. Simple components


Consists of 17 simple tools to solve all cases

03. Intuitive fixture removal protocol

Simple yet reliable performance, as it is completed with just a reverse rotation after attaching the Universal Driver.

04. Various screw removal protocols

Applications

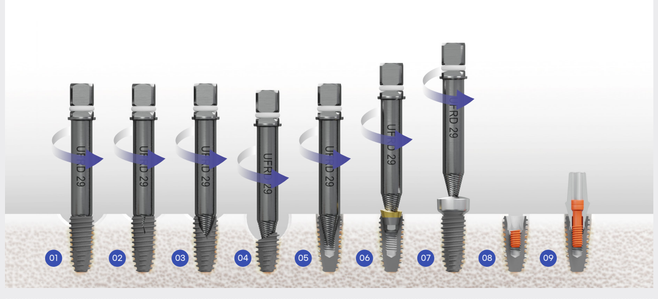

01 Adjacent Bone Loss
02 Crack
03 Partial Fracture (Vertical)
04 Complete Fracture (Horizontal)
05 Damage to the internal structure of the Fixture
06 Abutment Fracture
07 Wear of Healing/Cover Screw Hex
08 Abutment Screw Fracture
09 Wear of Abutment Screw Hex

D:O IMPLANT
www.dioimplant.com

New Dr. SOS+ Kit


One Kit Solution for Implant Failure Cases



D:O New Dr. SOS+

DIO New Dr. SOS+ Kit
Components

Kit Code I DRS 05

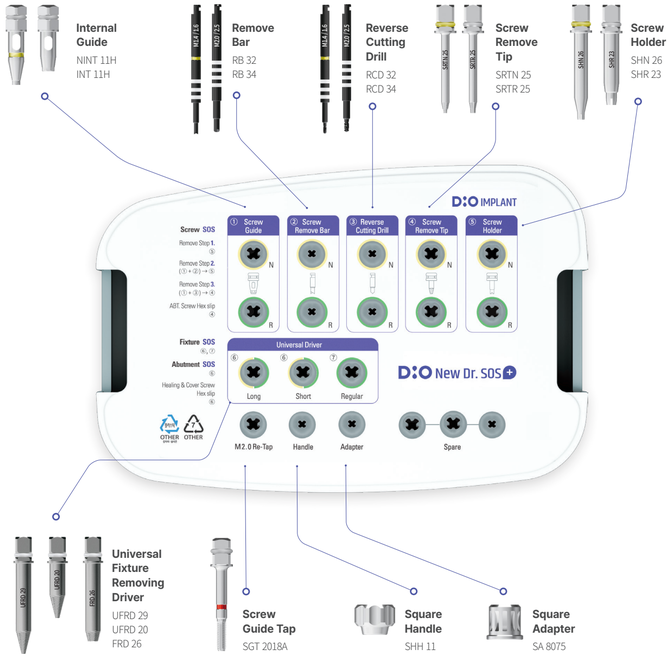

Square Wrench
SSW 050



*Bottom Panel Components

* SOS cases: Fractured Implants, Abutments, Screws and damaged hex slots, etc.

Rescue implant fracture

One tool to resolve all types of SOS cases!


01 Removal of Fractured Implant

* CCW : Counterclockwise, CW : Clockwise



Implant fractured — Rotate **Universal Driver** in CCW — Replace with same or similar size implant

02 Removal of Fractured Abutment

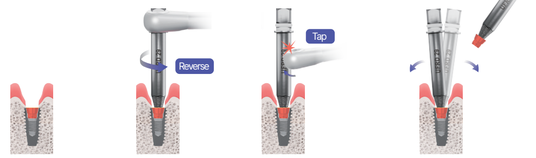

Abutment fractured — Rotate **Universal Driver** in CCW — Gently tap the handle with a **Torque Wrench** — Rock the **Universal Driver** sideways to remove abutment

03 Removal Healing Abutment / Cover Screw with Damaged Hex



Damaged Hex — Rotate **Universal Driver** in CCW — Remove Healing Abutment or Cover Screw


How to Use

Removal of Fractured Abutment Screw

Fractured Abutment Screw Can Be Removed With a High Successful Rate


04 Removal of Fractured Abutment Screw

First try Use Screw Holder

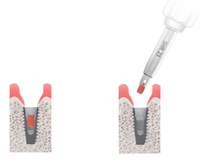

Fractured screw is clogged in fixture screw hole — Pick up with **Screw Holder** and remove
* Hand Mode

Second try When you fail in your first try, use the Screw Remover Bar and the Screw Holder

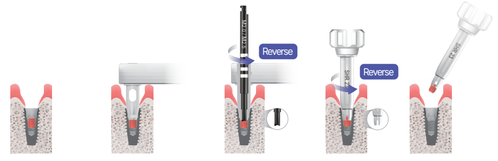

First try failed — Place Guide at Fixture Connection — Use Screw Remove Bar rotating in CCW
*Recommended speed : 20 rpm — Pick up with **Screw Holder** and remove
* Hand Mode


Third try When you fail in your second try, use the Screw Remover Bar and the Screw Holder



Second try failed — Place Guide at Fixture Connection
*Make sure connection is secure — Use Reverse Cutting Drill rotating in CCW
*Recommended speed: 1,200 rpm
*Must use with irrigation & pumping
*Drill at least 1mm deep — Place **Screw Remove Tip** with Guide connected, pick up screw and remove
* Hand Mode

How to Use

05 Removal of Abutment Screw with Damaged Hex

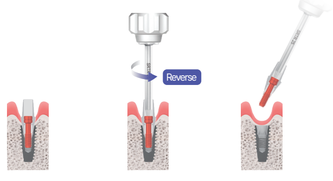

Damaged Hex — Use **Screw Remove Tip** rotating in CCW — Remove Abutment Screw

06 Restoration of damaged threads in fixture hole



Hole threads damaged — Use appropriate Re-Tap or fixture threads, turn clockwise
* Hand Mode
(Contact sales rep for selecting Re-Tap)

Reference Detach the Fixture from the Driver



Insert the fixture in the X-shaped hole on the **Square Wrench.**
Push down and turn clockwise to separate the Universal Driver.

How to Use